


Underivatized Cold Medications (Basic Drugs) Rtx®-35 Amine

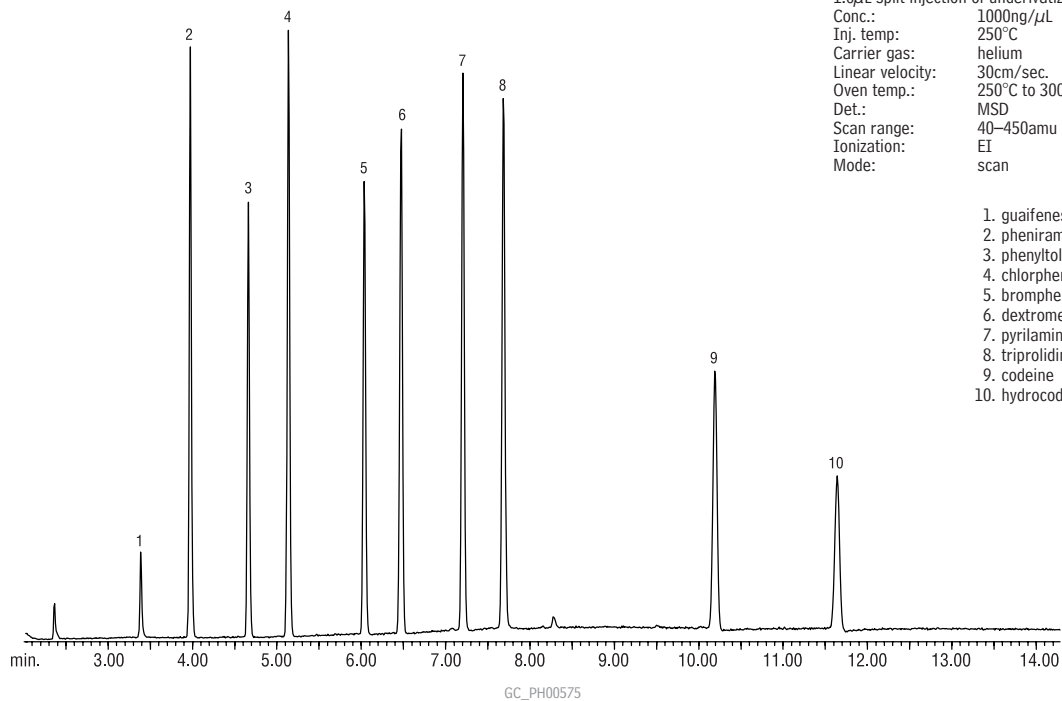

Rtx®-35 Amine 30m, 0.25mm ID, 0.50 μ m (cat.# 11338)
1.0 μ L split injection of underivatized cold medicines
Conc.: 1000ng/ μ L
Inj. temp: 250°C
Carrier gas: helium
Linear velocity: 30cm/sec.
Oven temp.: 250°C to 300°C @ 7°C/min. (hold 7 min.)
Det.: MSD
Scan range: 40–450amu
Ionization: EI
Mode: scan

1. guaifenesin
2. pheniramine
3. phenyltoloxamine
4. chlorpheniramine
5. brompheniramine
6. dextromethorphan
7. pyrilamine
8. triprolidine
9. codeine
10. hydrocodone